


PRODUCT

**URINE COLLECTION BAG SM7 "CLOSED SYSTEM" STERILE
 PROVIDE WITH DRIP CHAMBER**

FEATURES

- Bag's body:** made of atoxic, soft and strong PVC.
- Bag's colour:** dull white and transparent.
- Bag's tubing:** large sectioned, made of atoxic and transparent PVC; PVC granule's high quality gives a high "memory" effect to the tubing making it very folding enduring. The external diameter is 9.5 mm. The tubing length could be adapted according to customers' need (between 90 an 150 cm).
- Connector:** universal, white, made of atoxic material (ABS); provided with sampling port.
- Sampling port:** latex free; it allows to effect drawings of urine samples; "needless" sampling port, which guarantees optimal functionality for 10 (ten) drawings.
- Non-return valve:** close to prevent urine and gases backflow.
- Tap and emptying:** by the "cross-valve" tap, it allows the bags to be easily and quickly emptied with the guarantee of preventing obstructions during the outflow.
- Duration use:** 5 days. The "T" tap assures hygiene and l'use till 8-10 days if used in connection with the auxiliary emptying bag "Vià!".
- Bag's capacity:** 2 l.
- Accessories:** plastic hanger assemble on the bag;
 drip chamber;
 air filter for the compensation of the pressure;
 clamp to help the operator during samples drawing.
- Sterilization:** made by ETO and CO2

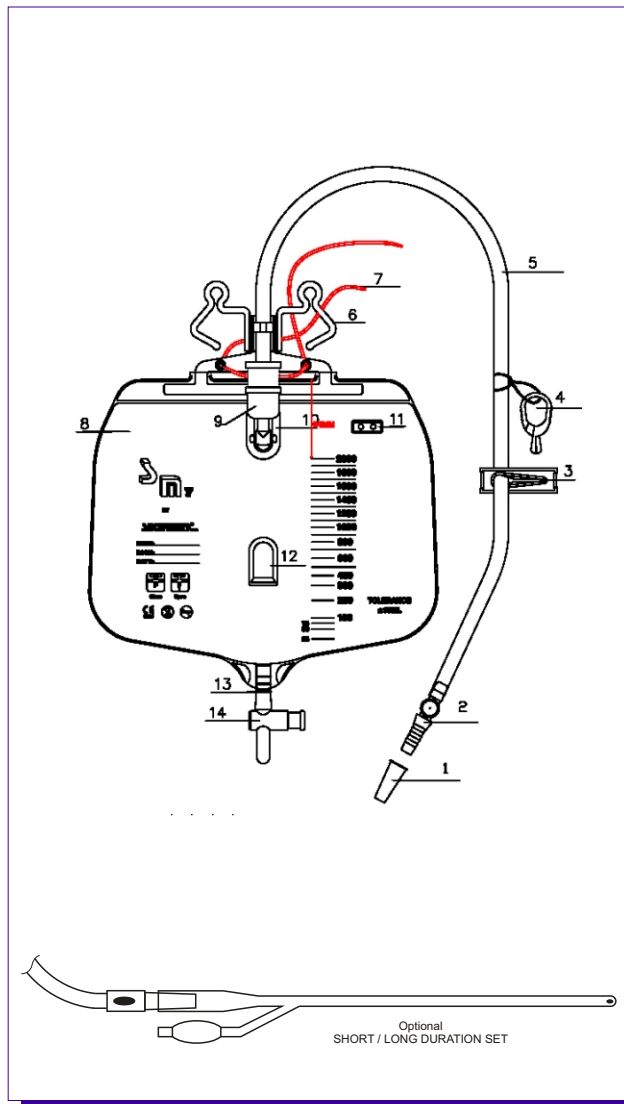

CODE	DESCRIPTION	PIECES PER CARTON	CARTON WEIGHT	CARTON DIMENSIONS	CARTON VOLUME
004011E17903	SM7 urine .bag, sterile, closed system, 2 l., tube of 90 cm and clamp, with drip chamber and air filter and sampling port connectors "needless", with folding lateral "cross-valve" tap and assemble plastic hanger.	20	5,6 Kg.	0.39 x 0.59 x 0.25 m.	0.057 m. ³
004011E17C23	SM7 urine .bag, sterile, closed system, 2 l., tube of 120 cm and clamp, with drip chamber and air filter and sampling port connectors "needless", with folding lateral "cross-valve" tap and assemble plastic hanger.	20	5,9 Kg.	0.39 x 0.59 x 0.25 m.	0.057 m. ³